

Figure S1. Standard curve and amplification plot of ISGs and GAPDH.

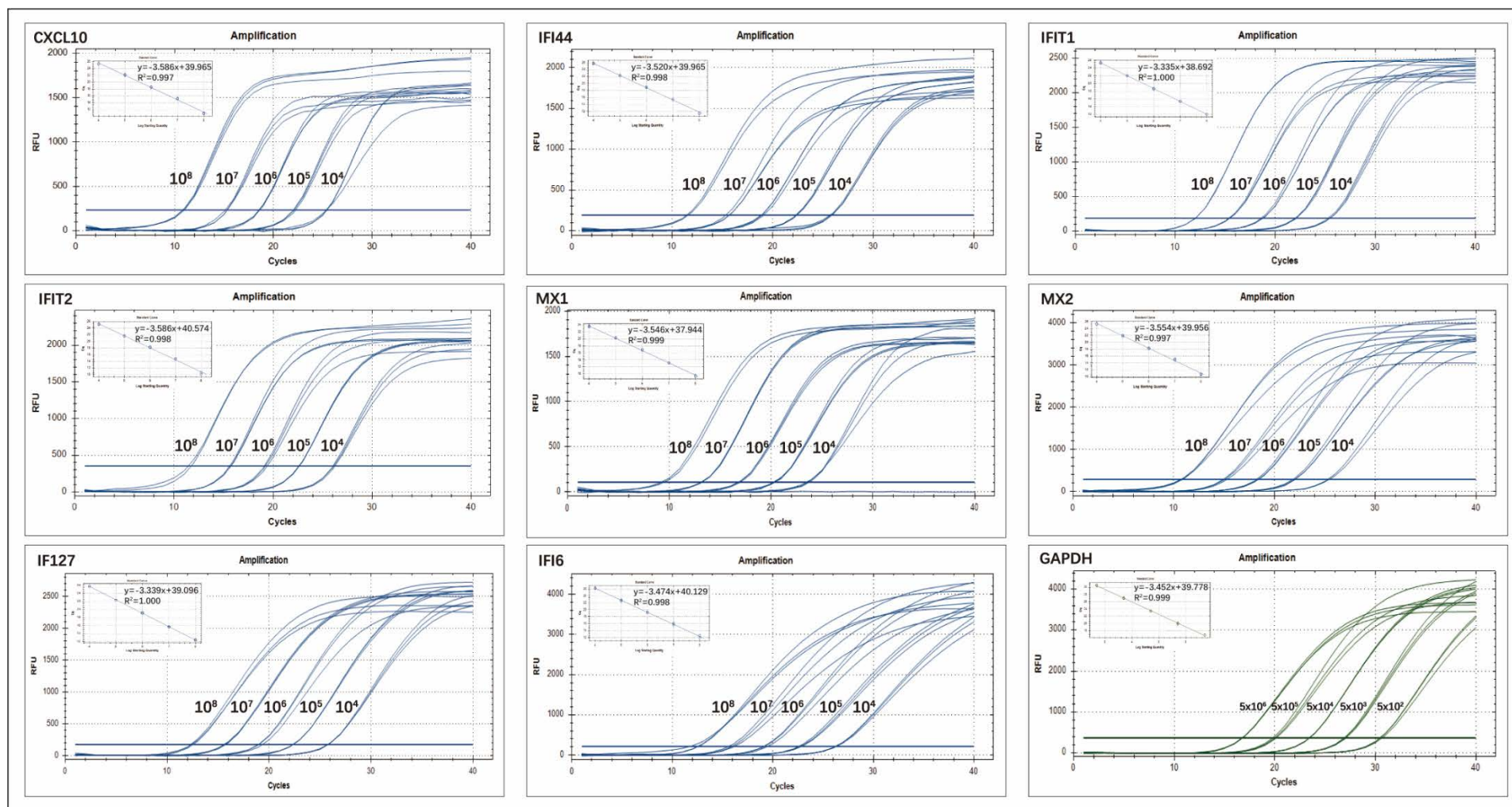


Table S1. Primers for amplifying ISG and GAPDH by RT-qPCR assay

Target ISG	Labeled	Sequence (5'→3')
CXCL10	Forward Primer	CTTTCTGACTCTAAGTGGCATTC
	Reverse Primer	CACCCTTCTTTTTCATTGTAGCAA
	Probe	FAM-CCTGCAAGCCAATTTTGTCCACGTGT-BHQ1
IFI44	Forward Primer	TGCCTTGAGAACTTATGAACCA
	Reverse Primer	ATATCCCAGTTGTATTAGTGCC
	Probe	FAM-AAAGCCTGATGCGTTACATGCCCTT-BHQ1
IFIT1	Forward Primer	TTGGTTTTAAGGAAACTTCGGAGA
	Reverse Primer	GCCCGCTCATAGTACTCCA
	Probe	FAM-AAAGCTTGAGCCTCCTTGGGTTTCG-BHQ1
IFIT2	Forward Primer	CAATGATAATCTCTTCCGTGT
	Reverse Primer	CAGTTGTTTCGCTACAGG
	Probe	FAM-ATTCTTGCCAGCCTCCATGCTC-BHQ1
MX1	Forward Primer	AGGCAAGGTCAGTTACCA
	Reverse Primer	CTCCAGGGTGATTAGCTC
	Probe	FAM-TTCCCATTCTTCCCCGGCGAT-BHQ1
MX2	Forward Primer	GGCAGAAAGACTTACCACT
	Reverse Primer	AACATCTTGTCGGCCTC
	Probe	FAM-CTGGTGGCTCTCCCTTATTTGTCCT-BHQ1
IFI27	Forward Primer	GGGAATCGCCTCGTCCTC
	Reverse Primer	GTAGAACCTCGCAATGACAGC
	Probe	FAM-CCCAGTGACTGCAGAGTAGCCACA-BHQ1
IFI6	Forward Primer	TGATGAGCTGGTCTGCGA
	Reverse Primer	ATCAGGGCACCAATATTACCTA
	Probe	FAM-CTGCCACCAGCCCCGAGG-BHQ1
GAPDH	Forward Primer	CAACAGCGACACCCAC
	Reverse Primer	ATGCCAGCCCCAGCGTCA
	Probe	HEX-AGCCAAATTCGTTGTCATACCAG-BHQ1